

Read Before Using AERIS Airway®

The Detachable Nasal Cannula is for pre-procedure and post-procedure oxygenation.





The AERIS Airway fits adult-size Guedel OPAs (80mm-110mm).

- AERIS® does not fit in Guedel OPAs with a nonstandard internal lumen, such as a catheter guide channel or other modified airway passage. Examples include the Teleflex® Hudson RCI® Cath-Guide® Guedel Airway and Intersurgical® Guedel Airway.

The AERIS Airway fits adult-size NPAs (24-36Fr). The following manufacturers' NPAs are recommended for the best fit:

- AirLife® SunMed PVC NPA (*remove the 15mm connector first*)
- Teleflex RUSCH Robertazzi NPA


The AERIS Airway® Hub Must Be Fully Seated in an OPA or NPA for Optimal O₂ Delivery and EtCO₂ Sampling

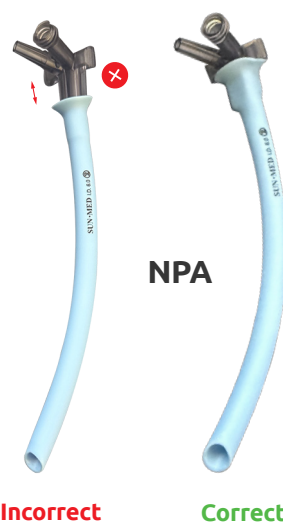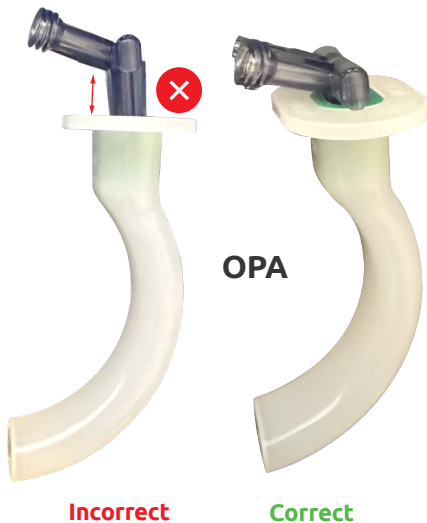

Top view (seated)




Clinical Essentials Before Using AERIS Airway®

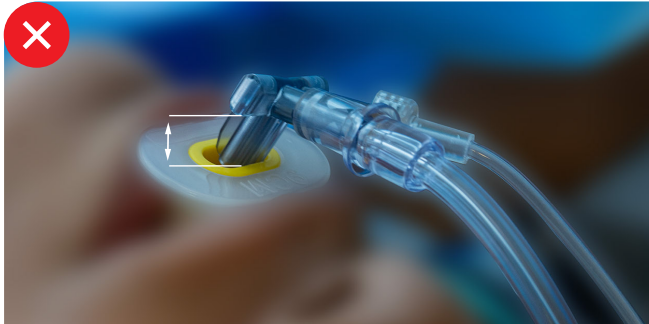


Incorrect

OPA

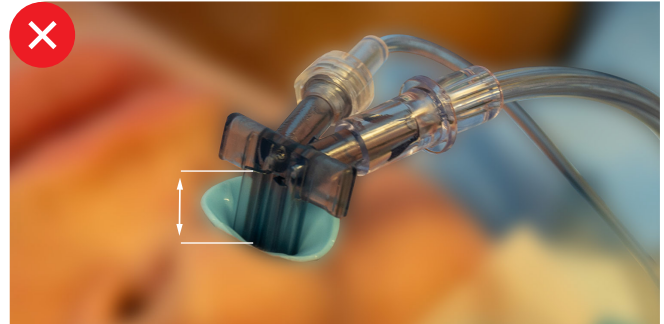


Incorrect

NPA

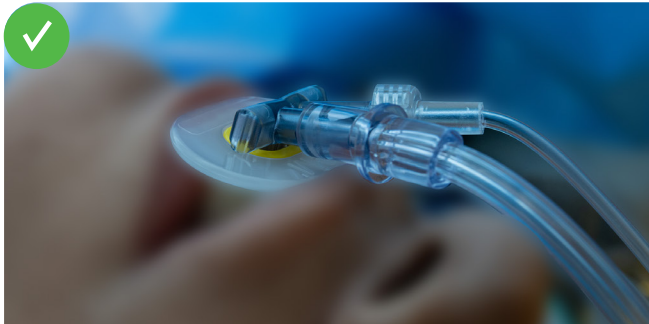

Correct

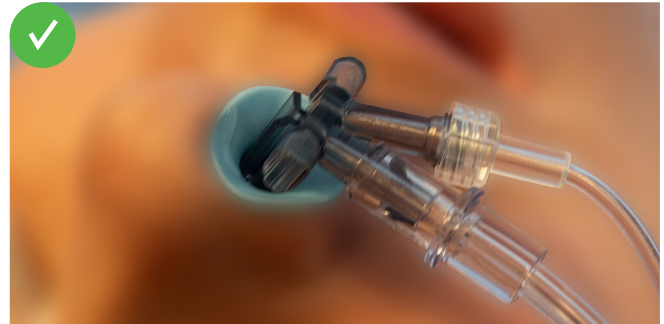

Correct

For More
New User Info:



844.772.3747
info@airwayinnovations.com

©2026 Airway Innovations